

CLEVELAND OPEN CUP FLASH POINT TESTER PTL2V-06



CLEVELAND OPEN CUP FLASH POINT TESTER


PTL2V-06

Petroleum testing is the analysis during upstream, midstream, and downstream production processes of petroleum products. It is most commonly used to test petroleum product, its product components, byproducts of crude oil, fuel, natural gas, upstream oil and gas and other formats of petroleum.

Used in Petroleum Industry, PVC Pipe Industry.

PTL2V-06 CLEVELAND OPEN CUP FLASH POINT TESTER

It adopts special heating furnace to ensure the safety of test. It is applicable to all petroleum products with flash points above 79°C and below 400°C except fuel oils. The heating power is continuously adjustable. It can meet requirements of test. The instrument is fully self-contained complete. The operator can do determination as long as connecting with coal gas or other civil gas. With fast cooling function, it can improve test efficiency.

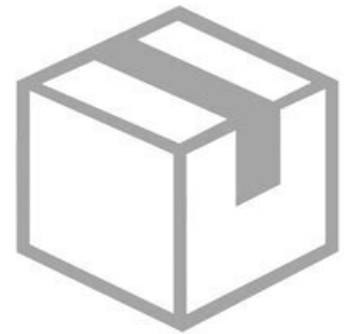

Product Image Coming Soon

SPECIFICATIONS

Model	PTL2V-06
Igniting device	
Ignition source	Coal gas (or civil gas)
Nozzle aperture	About 0.8 mm
Thermometer	(-6~400) °C. Scale is 2 °C
Ambient temperature	(15~35) °C
Relative humidity	≤ 85 %
Test flame applicator	It applies the test flame automatically
Heating device	furnace heating, no naked fire, explosion prevented
Power supply	AC (220±10%) V, 50Hz.
Power consumption	≤650W
Dimension	350×290×350 mm (thermometer is not included)
Net weight	5.5 kg

LABSTAC

Labstac LLC

82 Wendell Avenue, STE 100, Pittsfield, MA, 01201, USA
Email: contact@mail.labstac.com | Website: labstac.com