

# PARAFFIN TRIMMER PTR11-230



# PARAFFIN TRIMMER PTR11-230

Paraffin Trimmer is used to remove excess wax from the outside of a tissue cassette and to clear the edge of the cassette of excess wax after paraffin embedding. This device features a ridged heating surface where the temperature can be adjusted in a range from ambient to 90°C.

Used in Medical, Research, Tissue, Biopsy, Mega and Embedding.

Also known as Wax Trimmer, Laboratory Wax Remover, Laboratory Paraffin Remover.

## PTR11-230 PARAFFIN TRIMMER

Full computer automated control.


Utilizes new type of heating elements that enable fast heating, energy, efficiency and reliability.

Temperature sensors are made in USA, precise and reliable.

Automatically saves the latest operational settings.

Memory and automatic restoration function, the preset temperature can be automatically saved after device operation.

LED displays the heating status.

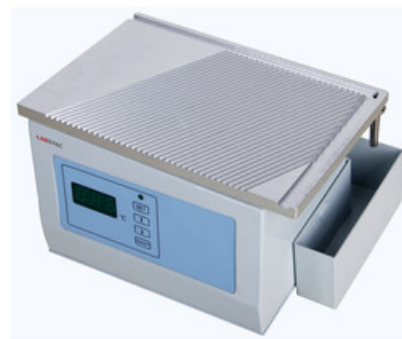

## SPECIFICATIONS

|                               |   |
|-------------------------------|---|
| Model                         | PTR11-230                               |
| Working Area                  | 230×160 mm                              |
| Temperature Range             | Continuously adjustable within 0 - 90°C |
| Temperature Control Precision | ±1°C                                    |
| Electrical Requirement        | 220V±10%50 Hz or 110V±10% 50/60 Hz      |
| Power                         | 150 W                                   |
| Dimension                     | 255W×180D×155H mm                       |
| Weight                        | 3 kg                                    |

# LABSTAC

**Labstac LLC**

82 Wendell Avenue, STE 100, Pittsfield, MA, 01201, USA

Email: [contact@labstac.com](mailto:contact@labstac.com) | Website: [labstac.com](http://labstac.com)